

Biotherm

Code: Biotherm


High protection thermic sheet against hypothermy
favours the recovery of the lost heat



Radiotransparent




Antistatic

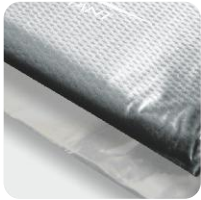

Non- allergic



Lear resistant

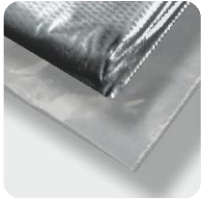

Non sliding



CE

biomedical



Biotherm

Code: Biotherm

High protection thermic sheet against hypothermy
favours the recovery of the lost heat

To order:

CODE	PACKAGE
B1000 Adults cap	120
B1010 Pediatric cap	120



B3000 Glove for pulsossimeter	120
-------------------------------	-----



B1100 Legging 100 cm	50
----------------------	----



B1150 Special pediatric sheet (75x77cm)	100
---	-----



B1200 Sheet S (60x60cm)	80
B1250 Special size sheet (60x90cm)	50
B1255 Special size sheet (90x120cm)	50
B1300 Sheet M (120x120cm)	50
B1400 Sheet L (210x120cm)	40



B0950 Roll (1,15x50m)	2
-----------------------	---



In the post-operative phase Biotherm facilitates the recovery of any heat lost, thereby easing the patient's awakening.

Biotherm consists of a support in polyethylene, treated so as to be antistatic, non-conducting and X-ray penetrable, and laminated at the polymer level with an aluminizing agent. Subsequently, it is embossed and heat-bonded, without the use of any solvents, with natural long-fiber cellulose non-woven pulp.

Biotherm can be sterilized, reflects heat but not light, and holds perspiration without eliminating it. It is non-slip, with good draping features, not irritating, non-allergic, odourless, tear-resistant and can be cut to any size without fraying or shedding fibres or particles.

It does not contain additives, it burns without emitting fumes, volatile substances or toxic and polluting residues, it wets on the inside in scattered spots and on the outside it is water-repellent. Biotherm is made in Italy.

M-721-521 rev. A del 03/09/2012



biomedical

Biomedical srl Via G.B.Lulli, 43 50144 Firenze Italy
Tel +39 055 352678 +39 055 332822 Fax +39 055 331629
www.biomedical-bio.net sales@biomedical-bio.net

Bone biopsy needles
Guillotine needles for biopsy on soft tissues
Cyto-aspiration needles
Different specialistic needles and sets
Loco-regional anaesthesia needles and sets
Thermic and absorbing sheets